

## Sacchrosomes - Human Cytochrome P450s in a Yeast Expression System.

# Human CYP2D6 + P450 Reductase

### *Product overview*

Catalogue Number CYP2D6-1      Lot Number 2D6-10-07


Human CYP2D6 is expressed in both hepatic and extra-hepatic tissues. It is known to be constitutively active and can metabolise 20-25% of all clinically used drugs.

Membranes consist of Human Cytochrome P450 2D6 and P450 Reductase enzymes bound within a yeast microsomal fraction.

### *Product Data*

Pack Size	0.5 nmol
Volume per tube	0.5 ml
Cytochrome P450 Content	119 pmol.mg <sup>-1</sup>
Protein Concentration	8.4 mg.ml <sup>-1</sup>
Specific Activity	0.187 pmol CHC.min <sup>-1</sup> . pmol P450 <sup>-1</sup>
Cytochrome P450 Reductase Activity	610.7 nmol MTT reduced.min <sup>-1</sup> .mg protein <sup>-1</sup>

**Fluorometric assay of CYP2D6- Graph depicts product formation (3-cyano-7 hydroxycoumarin) in pmol over time in min.**

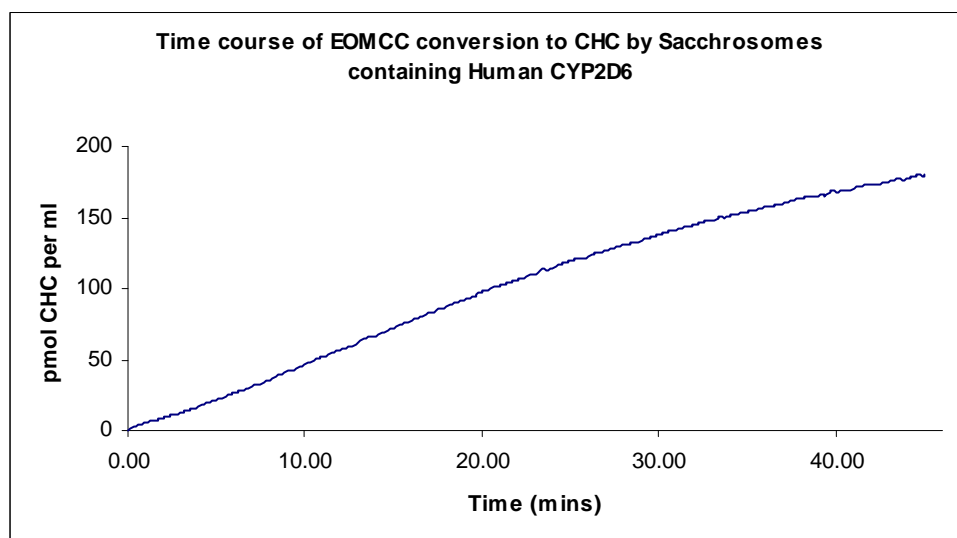

### **Specific Activity**

CYP2D6 activity assay performed in a microplate-based fluorometric assay with EOMCC (Invitrogen) as substrate in 0.1M phosphate buffer (pH 7.4). Excitation at 400nm and Emission at 460nm Temperature held at 37°C.

0.1 ml of reaction mix contains 1.3mM NADP<sup>+</sup>, 3.3 mM Glucose-6-phosphate, 3.3 mM Magnesium Chloride, 0.04U of Glucose-6-phosphate dehydrogenase and 10 µM of EOMCC. 2.5 pM of Sacchrosome 2D6 is added per reaction.

The conversion of EOMCC substrate to 3-cyano-7-hydroxycoumarin (CHC) product is measured over time. Values converted using a standard curve of 3-cyano-7-hydroxycoumarin.

### **Cytochrome P450 Content**

CO binding assay performed in a cuvette format using a dual beam spectrophotometer scanning from 500 to 400nm. Spectral difference of microsomes measured in a phosphate glycerol buffer with the addition of sodium dithionite with and without CO perfusion.

### **Cytochrome P450 Reductase Activity**

Reduction of MTT by Cytochrome P450 Reductase utilising a regenerating system in a phosphate buffer was measured over time.

### **Protein Concentration**

Total protein was measured using a microplate-based Bradford assay method with BSA as a standard.

### **Product Use**

For best stability thaw on ice, aliquot suitable quantities for your studies and store at -80°C.

Microsomes are supplied in a buffer containing Water, Tris, EDTA and Glycerol which are unlikely to interfere with most assays.

Studies indicate product stability at -80°C for at least 12 months.

### **Safety**

This product is not suspected to contain any pathogenic or hazardous materials. However, since these properties have not been investigated handle with care in accordance to your normal laboratory practices.

This product is only intended for *in vitro* research use and is not licensed as a drug, therapeutic or diagnostic tool for humans or animals.

Sacchrosomes™ of DMBioscience Ltd, Patents pending, ©.